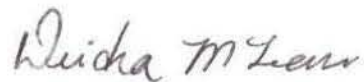


The author hereby certifies that the use of any copyrighted material in the thesis manuscript entitled:

**"Enamel Bond Strength of New Universal Adhesive Bonding Agents"**

is appropriately acknowledged and, beyond brief excerpts, is with the permission of the copyright owner.

A handwritten signature in dark ink, appearing to read "Deidra McLean". The script is cursive and fluid.

**Major Deidra McLean**  
**Advanced Education in General**  
**Dentistry Residency**  
Uniformed Services University  
Date: **6 June 2013**

**Uniformed Services University  
of the Health Sciences**

**Manuscript/Presentation Approval or Clearance**

**INITIATOR**

1. USU Principal Author/Presenter: Maj Deidre McLean
2. Academic Title: Resident, Advanced Education in General Dentistry Residency (AEGD-2)
3. School/Department/Center: Air Force Postgraduate Dental School (AFPDS)
4. Phone: 210-671-9822
5. Type of clearance: ☒ Paper ☐ Article ☐ Book ☐ Poster ☐ Presentation ☐ Other
6. Title: Enamel Bond Strength of New Universal Bonding Agents
7. Intended publication/meeting: *Journal of the Canadian Dental Association*
8. "Required by" date: 15 July 2013
9. Date of submission for USU approval: 7 June 2013

**CHAIR OR DEPARTMENT HEAD APPROVAL**

1. Name: Kraig S. Vandewalle, Director, Dental Research
2. School/Dept.: AFPDS / AEGD-2
3. Date: 6 June 2013

**\*Note:** *It is DoD policy that clearance of information or material shall be granted if classified areas are not jeopardized, and the author accurately portrays official policy, even if the author takes issue with that policy. Material officially representing the view or position of the University, DoD, or the Government is subject to editing or modification by the appropriate approving authority*

Chair/Department Head Approval: Kraig A. Vandewalle Date 6 June 2013

---

### SERVICE DEAN APPROVAL

1. Name: Col Thomas Schneid



2. School (if applicable): AFPDS

3. Date:

4. ☐ Higher approval clearance required (for University-, DoD- or US Gov't-level policy, communications systems or weapons issues review").

**\*Note:** *It is DoD policy that clearance of information or material shall be granted if classified areas are not jeopardized, and the author accurately portrays official policy, even if the author takes issue with that policy. Material officially representing the view or position of the University, DoD, or the Government is subject to editing or modification by the appropriate approving authority.*

### COMMANDER APPROVAL

1. Name:



2. School (if applicable):

3. Date:

4. ☐ Higher approval clearance required (for University-, DoD- or US Gov't-level policy, communications systems or weapons issues review").

**\*Note:** *It is DoD policy that clearance of information or material shall be granted if classified areas are not jeopardized, and the author accurately portrays official policy, even if the author takes issue with that policy. Material officially representing the view or position of the University, DoD, or the Government is subject to editing or modification by the appropriate approving authority.*

### PDC DEAN APPROVAL

1. Name:

2. School (if applicable):

3. Date:

4. \_\_Higher approval clearance required (for University-, DoD- or US Gov't-level policy, communications systems or weapons issues review").

**\*Note:** *It is DoD policy that clearance of information or material shall be granted if classified areas are not jeopardized, and the author accurately portrays official policy, even if the author takes issue with that policy. Material officially representing the view or position of the University, DoD, or the Government is subject to editing or modification by the appropriate approving authority.*

---

Dean/VP Signature/Date

#### **VICE PRESIDENT FOR EXTERNAL AFFAIRS ACTION**

1. Name:

2. Date:

3. \_\_USU Approved or

\_\_DoD Approval/Clearance required

4. \_\_Submitted to DoD (Health Affairs) on (date):

Or \_\_Submitted to DoD (Public Affairs) on (date):

5. \_\_DoD approved/cleared (as written) or \_\_DoD approved/cleared (with changes)

6. DoD clearance/date:

7. DoD Disapproval/date:

---

External Affairs Approval

Date

## **Enamel Bond Strength of New Universal Bonding Agents**

### **Abstract**

The demand for esthetic, conservative restorations has inspired considerable research and marketing by manufacturers to develop effective dental adhesives that are faster and easier to place and less technique sensitive. In an attempt to simplify the bonding procedure, the bond strength of some of the simplified materials has been shown to be reduced. Concerns about hydrolytic degradation of the bonded interface have also been raised. Most recently, two new Universal Bonding Agents have been marketed for use as self-etch or etch-and-rinse adhesives depending on the dental substrate and clinician's preference. Objective: The purpose of this study was to evaluate the shear bond strength of composite to enamel using two new universal adhesives compared to a self-etch adhesive when applied in self-etch and etch-and-rinse modes after 24 hours and 6 months storage. Materials and Methods: One-hundred twenty enamel sections were obtained from human third molars and mounted in PVC pipe. A small area of the enamel was flattened using a diamond wheel and 600-grit silicon-carbide paper. The enamel specimens were randomly divided into 3 groups of forty and assigned to two Universal adhesives, All-Bond Universal (Bisco) and Scotchbond Universal (3M/ESPE), and one self-etch adhesive, Clearfil SE (Kuraray). The Universal and self-etch adhesive groups had twenty specimens bonded in a self-etch mode and twenty specimens bonded in an etch-and-rinse mode with a separate 34% phosphoric acid etch. Each enamel specimen had its assigned adhesive applied as per manufacturers' instructions, then placed in an Ultradent Jig. Composite (Filtek Z250, 3M/ESPE) was inserted into the mold in three increments to a height of 4 mm and cured for 20 seconds per increment. The 6 groups were then divided into two subgroups with ten specimens each. One subgroup was stored for 24 hours and the second subgroup for 6 months in 37°C distilled water and tested in shear in a universal testing

machine (Instron). A mean shear bond strength value (MPa) and standard deviation were determined per group. Data was analyzed with 3-way ANOVA/Tukeys ( $\alpha=0.05$ ). Following testing, each specimen was also examined using a 10X stereomicroscope to determine failure mode as either: fracture at the adhesive interface, cohesive fracture in composite, enamel, or dentin, or mixed fracture involving a combination of adhesive and cohesive fractures. Results: The 3-way ANOVA found a significant difference between groups based on bonding agent ( $p<0.001$ ) and surface treatment ( $p<0.001$ ) but not on time ( $p=0.943$ ) with no significant interaction ( $p>0.05$ ). Clearfil SE in etch-and-rinse and self-etch modes had more mixed fractures than either All Bond Universal or Scotchbond Universal in either mode. Conclusions: Etching enamel significantly increased the shear bond strength of composite to enamel. Clearfil SE had significantly greater bond strength to enamel than Scotchbond Universal and All-Bond Universal, which were not significantly different from each other.

## **Introduction**

Adhesive dentistry has been around for over fifty years since it was first introduced by Buoncore in 1955.<sup>1</sup> Since that time, there has been a constant evolution in the field of adhesive dentistry with the progressive introduction of seven generations of adhesive bonding agents.

Adhesive bonding agents must be capable of providing equally effective bonds to both enamel and dentin despite being vastly different structures in terms of composition and natural variability. Enamel's composition is primarily inorganic (86%) hydroxyapatite with 2% organic content and 12% water while dentin is composed of 50% inorganic mineral, 30% organic collagen and 20% water.<sup>2,3</sup> Enamel is a homogenous structure, while dentin is highly variable depending on several factors including age, dentinal

tubule number and size, and previous exposure to carious, chemical or mechanical stimuli.<sup>2</sup>

Other variables that may interfere with adhesive bonding include the creation and removal of a smear layer, as well as its thickness. Dentinal wetness may also affect bond strength if the tooth is left too wet or too dry following acid etching.<sup>4</sup> Matrix metalloproteinases also affect adhesive bonding over time.<sup>5</sup> Other challenges to adhesive dentistry in addition to differences between enamel and dentin include the presence of moisture in the working area, technique sensitivity of the materials, biocompatibility of materials, the requirement for a gap free restorative interface, and the requirement for the bonding agents to rapidly develop high bond strengths.

The basic mechanism of adhesion between tooth structure and composite resin is based on an exchange process. Minerals from hard tissue are replaced by resin monomers that effectively create a micromechanical bond.<sup>6</sup> Despite the similarities between adhesives, the composition of these materials and the manner in which they are applied differ. The demand for simpler, more user-friendly and less technique-sensitive adhesives has inspired manufacturers to develop new adhesives at a rapid rate.<sup>7</sup>

Currently, there are four generations of dental adhesives available to dentists including 4th, 5th, 6th and 7th generation adhesive bonding agents. In addition to the generation classification, there is also an adhesive classification system. This hierarchy classification system includes two major categories of adhesives: etch-and-rinse adhesives and self-etch adhesives. These major categories are further divided into four subtypes: three-step etch-and-rinse, two-step etch-and-rinse, two-step self-etch and one-step self-etch. The two-step etch-and-rinse and one-step self-etch are also referred to as simplified adhesives because the primer and adhesive are combined. The one-step self-etch adhesives may be further subdivided into 'two-component' and 'single-

component' one-step adhesives.<sup>7</sup> These classification systems and how they relate are demonstrated in Fig 1.

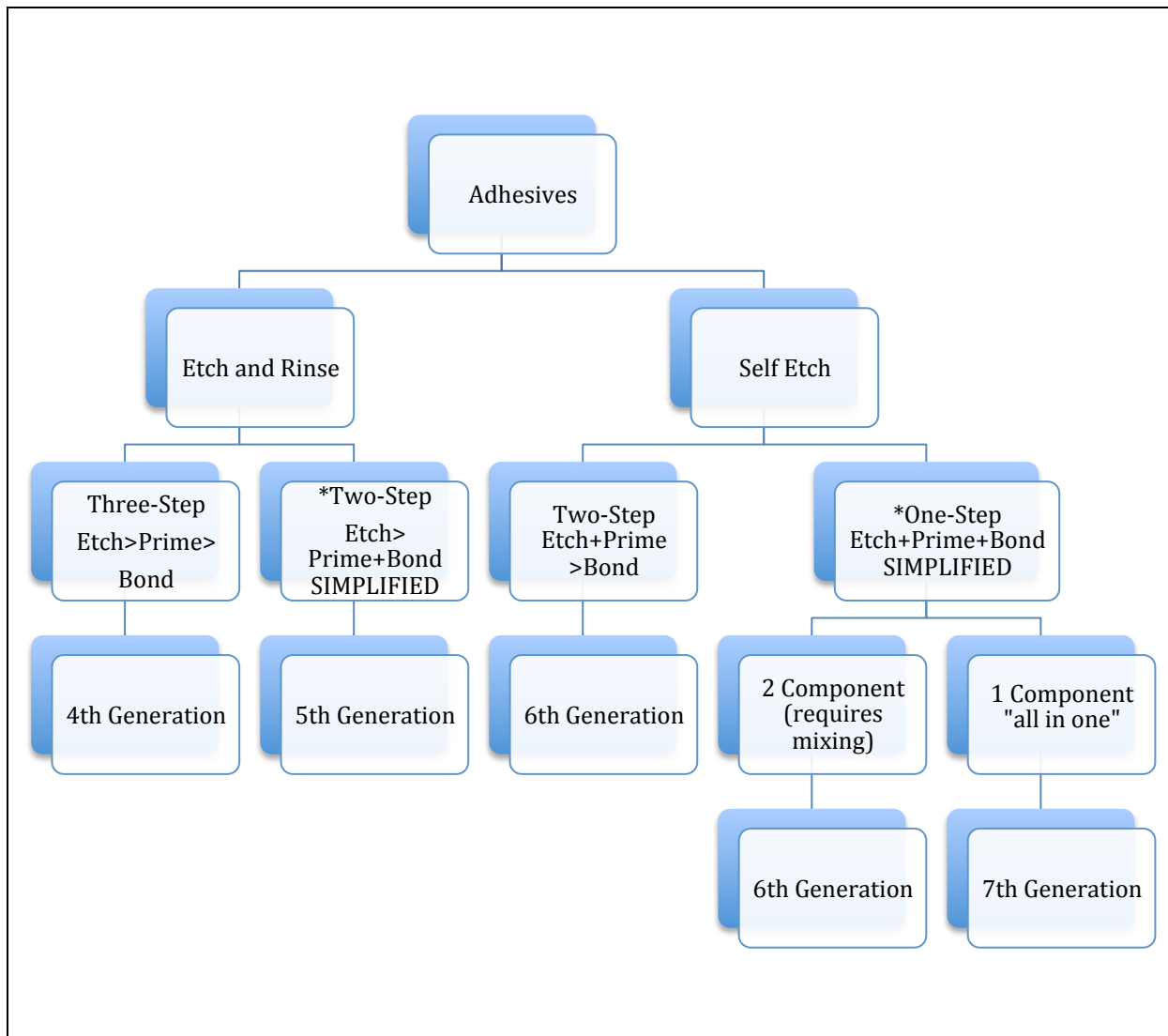


Fig 1. Classification of Adhesive Bonding Agents

Fourth generation or three-step etch-and-rinse adhesive bonding agents were developed in the early 1990's and are considered multi-step adhesives involving three separate applications including acid etching, application of the primer, followed by application of a separate adhesive. Fifth generation or two-step etch-and-rinse or simplified etch-and-rinse adhesives involve acid etching, followed by the combined application of a primer and an adhesive. The 6th generation or two-step self-etch



adhesives involve application of an acidified primer followed by application of the adhesive resin. The one-step self-etch adhesives, also known as the simplified self-etch adhesives, involve application of a combined acidified primer, and the adhesive resin in a single step. The two-component one-step self-etch adhesives, which are also 6th generation adhesive bonding agents, separate the active ingredients. Specifically, the functional monomer is separated from water, theoretically providing a longer shelf life, but additional and adequate mixing of both components is required. The single-component one-step adhesives, also known as 7th generation adhesive bonding agents, can be considered as the only true 'all in one' adhesives, combining the acidified primer and the adhesive resin and do not require mixing prior to application.<sup>7</sup>

Despite the various generations or adhesive classifications, there are significant differences between adhesive bonding agents even within the same class. For example, self-etch adhesives may vary greatly in their level of acidity. They may have strong, intermediately strong, mild or ultra-mild acid etchants.<sup>8</sup> Therefore clinical performance is highly product-dependent.

From the literature, Heintze<sup>9</sup> conducted the Meta analysis in 2010 that looked at the retention rates of cervical composite restorations bonded with various adhesive-bonding agents. As a result of these numerous clinical studies, it was concluded that the highest retention rates were achieved with the two-step, self-etch adhesive, Clearfil SE Bond (Kuraray, New York, NY), followed closely by the three-step etch-and-rinse adhesive, Optibond FL (Kerr, Orange, CA). Clearfil SE Bond had been shown to produce lower bond strength to enamel, particularly uncut enamel, but by selectively etching the enamel with phosphoric acid, Peumans<sup>10</sup> demonstrated that retention rates of Class V restorations after five years were 100% as opposed to 98% retentive without the selective etch step, but not statistically different. Studies have also shown significantly less marginal defects and staining with selective etching of enamel.<sup>10, 11</sup> Van

Meerbeek<sup>7</sup> also advocates the use of the selective etch technique using phosphoric acid on enamel because it produces the most durable bond to enamel that effectively seals and protects the more vulnerable bond to dentin against degradation.

A study by Peumans<sup>12</sup> that looked at the average annual failure rate of Class V composite resin restorations bonded with various dental adhesives revealed the three-step etch-and-rinse and two-step self-etch bonding agents to be most effective with a 4.8% and 4.7% annual failure rate respectively. The simplified adhesives, including the two-step etch-and-rinse and one-step self-etch adhesives had the highest annual failure rates of 6.2% and 8.1% respectively. According to a study by De Munck<sup>13</sup>, after approximately three months, all categories of dental adhesives start to exhibit mechanical and morphological evidence of bond degradation. The three-step etch-and-rinse adhesives were said to remain the 'gold standard' in terms of bond durability followed closely by the two-step self-etch adhesives. Any kind of simplification in the clinical application procedure resulted in loss of bonding effectiveness due to hydrolysis and elution of interface components.<sup>13</sup>

In late 2011 and early 2012, 3M ESPE and Bisco introduced two new Universal Bonding Agents. According to the manufacturers, these products can be used as etch-and-rinse, self-etch, and selective-etch adhesives for bonding to enamel or dentin as well as many indirect restorative surfaces depending on the clinician's preference. Reportedly, neither product requires refrigeration and can be stored at room temperature for two years.

The purpose of this *in vitro* study was to examine the shear bond strength of the new Universal Bonding Agents over time to enamel surfaces when used as an etch-and-rinse and self-etch adhesive compared to a two-step self-etch adhesive used in similar modes. The null hypothesis to be tested was that there would be no significant

difference in the shear bond strength of composite to enamel based on type of bonding agent, type of surface treatment, or time.

## **Materials and Methods**

Extracted human permanent third molars were stored in 0.5% Chloramine T solution for up to six months before being utilized. The teeth were visually examined and discarded if the enamel had caries or visible fracture lines. The crowns of the teeth were sectioned mesio-distally, then buccal and lingual sections were obtained by sectioning the crowns at the cemento-enamel junction using a water-cooled diamond saw (Isomet 5000, Buehler, Lake Bluff, IL). Each enamel specimen was mounted in polyvinylchloride (PVC) pipe using dental stone and bisacryl resin. After the stone had set, a small area of the enamel was cut flat using a diamond wheel bur then smoothed using 600-grit silicon-carbide paper.

The enamel specimens were divided into twelve groups with ten specimens each in order to compare the shear bond strength of different adhesives over time as depicted in Table 1. The adhesives that were compared included Clearfil SE (Kuraray), applied as a two-step self-etch and as a three-step etch-and-rinse adhesive; Scotchbond Universal Adhesive (3M ESPE) as a one-step self-etch adhesive and as a two-step etch-and-rinse adhesive and All Bond Universal as a one-step self-etch adhesive and as a two-step etch-and-rinse adhesive. For the adhesives applied with an etch-and-rinse technique, 34% phosphoric-acid gel etchant (Kerr Dental) was applied to the cut enamel for 15 seconds, rinsed with water for 15 seconds, then lightly air dried for three seconds before the application of the adhesive to the flattened enamel specimens as per manufacturer's instructions. The adhesives applied with a self-etch technique were applied directly to the cut enamel surfaces as per manufacturer's instructions. All adhesives were light cured with a light-curing unit (Bluephase 16i, Ivoclar Vivadent,

Amherst, NY) for 20 seconds. Irradiance was determined with a radiometer (LED Radiometer, Kerr Dental) and was considered acceptable if greater than 1200 mW/cm<sup>2</sup>.

**Table 1 Adhesive agents, surface treatments and storage times**

<b>Dental Adhesive</b>	<b>Immediate Group (24 hours)</b>	<b>Aged Group (6 month)</b>
Clearfil SE (self-etch)	CF SE 24 hr	CF SE 6 mo
Clearfil SE (etch-and-rinse)	CF E&R 24 hr	CF E&R 6 mo
Scotchbond Universal (self-etch)	SB SE 24 hr	SB SE 6 mo
Scotchbond Universal (etch-and-rinse)	SB E&R 24 hr	SB E&R 6 mo
All Bond Universal (self-etch)	AB SE 24 hr	AB SE 6 mo
All Bond Universal (etch-and-rinse)	AB E&R 24 hr	AU E&R 6 mo

Following application of the adhesives, the bonded specimens were placed in a jig (Ultradent Products, South Jordan, UT) and secured beneath a white plastic mold. The bonded area was limited to the 2.4 mm circle determined by the mold. Z250 (Kerr) composite resin was applied in three increments to a height of 4 mm. Each increment was polymerized for 20 seconds as recommended by the manufacturer using the light-curing unit. The immediate and aged shear bond strength specimens were stored for 24 hours and six months respectively in distilled water at 37°C in a laboratory oven (Model 20GC, Quincy Lab, Chicago, IL).

The shear bond strength of the specimens was tested in shear mode with a customized probe (Ultradent Products) in a universal testing machine (Model 5943, Instron, Norwood, MA) using a crosshead speed of 1.0 mm/min until failure. Shear bond strength in megapascals (MPa) was calculated from the peak load of failure in Newtons divided by the specimen surface area. The mean and standard deviation were determined per group. Data was analyzed with a three-way ANOVA and Tukey's post-hoc test ( $\alpha = 0.05$ ) to evaluate the effects of bonding agent (3-levels), surface treatment (2-levels), and time (2-levels) on the shear bond strength of composite to enamel.

Following testing, the specimens were examined under the microscope at 10X magnification to determine the failure mode as either: 1) adhesive fracture at the adhesive interface, 2) cohesive fracture in the composite, enamel, or dentin, or 3) mixed fracture involving a combined adhesive and cohesive fracture.

## **Results**

A 3-way ANOVA and Tukeys post-hoc tests revealed significant differences in the mean shear bond strengths of adhesive agents ( $p < 0.001$ ) and surface treatment ( $p < 0.001$ ) but not on time ( $p > 0.05$ ) and with no significant interaction ( $p = 0.943$ ).

In general, the bond strength of composite to enamel was significantly greater using Clearfil SE compared to Scotchbond Universal or All-Bond Universal, which were not significantly different from each other. Etching the enamel significantly improved bond strengths of the Universal adhesives compared to self-etching only. Storage time did not significantly affect bond strengths (Fig 2). A high percentage of mixed fractures including dentin corresponded to the higher bond strength values as found with Clearfil SE. The lowest bond strengths and the most adhesive failures occurred with All Bond Universal followed by Scotchbond Universal in self-etch mode. More mixed fractures were found for both Universal adhesives when applied in an etch-and-rinse mode. Storage time did not affect fracture mode (Fig 3).

## **Discussion**

This *in vitro* study demonstrated that the etch-and-rinse or selective-etch technique is an effective approach to achieving more predictable and stable micro-mechanical bonding of composite to enamel. However, this study also demonstrated that there is considerable variation between dental bonding agents and ultimately, the shear bond strengths produced, are largely material dependent.

Surface treatment significantly affected the shear bond strength of composite to enamel for the Universal bonding agents therefore the null hypothesis that there would be no difference based on surface treatment must be rejected.

From the results of this study, the shear bond strengths of the Universal adhesives to enamel were improved when the bonding agents were applied as two-step etch-and-rinse adhesives rather than one-step self-etch adhesives. This was attributed to an improved micromechanical bond being produced with the addition of the etch-and-rinse or selective-etch surface treatment. Etch-and-rinse or selective-etch adhesive systems are characterized by an initial etching step, typically with 32-37% phosphoric acid, followed by a thorough rinsing procedure that is responsible for the complete removal of the smear layer and selective dissolution of the enamel rods. This creates microporosities in the enamel that are readily penetrated by bonding agents via capillary attraction.<sup>14</sup> Following polymerization, micromechanical interlocking of tiny resin tags within the etched enamel surface provide a strong micro-mechanical bond to enamel.<sup>15</sup> The alternative self-etch approach only dissolves the smear layer but does not remove it, as there is no rinsing step, leaving the dissolved products to become incorporated within the bonded layer.<sup>16</sup> Furthermore, the degree of demineralization produced by self-etch adhesives depends largely on the acidity or etching aggressiveness of the functional monomer and is material dependent. According to Sunfield et al., the penetration of the adhesive system may be restricted to the more superficial enamel layers with creation of shorter resin tags when self-etch adhesives are used without a selective-etch step.<sup>17</sup> Erickson also found improved bond strengths with a pre-etch step and attributed this to the degree of etching or the etch morphology achieved.<sup>18</sup> When used without a selective or pre-etch step, even the most acidic of the self-etch adhesives only produced an etch pattern primarily involving the ends of enamel prisms with little effect on the interprismatic regions. The subsequent resin penetration was described as a negative

replica of the etch pattern with resin penetrating the etched prisms but not into the interprismatic unetched regions. The weakest acidic self-etch adhesives only achieved a fine pitting of the enamel surface and corresponding fine resin projections. Tay et. al also reported differences in the thickness of the enamel hybrid layers depending on the acidity of the adhesive and the resultant aggressiveness of apatite dissolution.<sup>19</sup>

Both Scotchbond Universal (pH=2.7) and All Bond Universal (pH=3.2) are considered ultra mild to mild acidic adhesives, therefore, the additional selective-etch step followed by thorough rinsing logically produced improved micromechanical bonds between the composite resin and the highly mineralized enamel substrate than without the selective-etch step. Nonetheless, neither the acidity of the adhesive agent, thickness of the hybrid layer, nor the length of the resin tags are solely responsible for bonding effectiveness and stability for all adhesives. This study confirmed previous studies and demonstrated that an ultra-mild (pH~2.7) self-etch adhesive, Clearfil SE Bond (Kuraray), was capable of achieving strong bonds to enamel with or without a selective-etch step.<sup>12,20</sup> This was particularly evident for the 6-month Clearfil SE groups in which the self-etch group produced the same mean shear bond strength as the etch-and-rinse group. The bonding effectiveness of Clearfil SE is believed to be related to the separation of the acidic monomers in its functional primer from its adhesive agent as well as its specific composition that includes methacryloxydecyl phosphate (MDP). The monomer, MDP contains phosphate groups, capable of producing ionic or chemical bonds with calcium in hydroxyapatite. The Universal Adhesives are ethanol and water based adhesive bonding agents and also contain MDP, however, by virtue of the etch, primer and adhesive components being combined, the bond strength may ultimately be reduced. The dihydrogenphosphate group in MDP can dissociate in water to produce two protons rendering it hydrolytically unstable.<sup>21</sup> Ultimately, the shear bond strength of a dental bonding agent is material dependent.

Within the limitations of this study, the bond strengths obtained confirmed significant differences in reliability of the adhesive bonding agents tested regardless of storage time and surface treatment. These differences are likely due to the specific chemical composition and formulation of each adhesive bonding agent. Therefore, the null hypothesis that there would be no significant difference in the shear bond strength of composite to enamel based on type of adhesive bonding agent must therefore be rejected. The Universal Bonding Agents are considered simplified adhesives and specifically as 5<sup>th</sup> or 7<sup>th</sup> generation bonding agents depending on their use with or without a selective etch step. As stated previously, one-step self-etch adhesives combine the acidified primer and adhesive agents and two-step etch-and-rinse adhesives combine the primer and adhesive and traditionally both have been more acidic and hydrophilic than the two-step self-etch adhesives that separate their acidic primers from the bonding agents. The hydrophilicity of the one-step self-etch adhesives has been stated to be the main disadvantage of these materials. This hydrophilicity leads to decreased bond strengths due to permeability of the adhesive layer and contributes to the hydrolysis of resin polymers and the consequent degradation of tooth-resin bonds over time.<sup>15, 22, 23</sup>

In terms of failure mode, Al-Salehi and Burke<sup>24</sup> reported that there is a relationship between the bond strength and fracture failure mode. From the results of this study, the higher bond strengths did correlate with greater mixed fractures or cohesive plus adhesive failure modes. Clearfil SE in etch-and-rinse and self-etch modes had more mixed fractures than either All Bond Universal or Scotchbond Universal in either mode. The Universal bonding agents produced more mixed fractures when used in etch-and-rinse mode than self-etch mode, which also correlated with bond strength. Storage time had no effect on failure mode.



From the results of this study, we failed to reject the null hypothesis that there would be no significant difference in the shear bond strength of composite to enamel based on time. Although the bond strengths of these new Universal adhesives were found to be inferior, the bond strengths of the materials between 24 hours and 6 months of water storage were not significantly different therefore longer storage times would be needed to determine the effect of bond strength over time.

## **Conclusions**

The new Universal bonding agents demonstrated higher shear bond strengths to enamel with the added selective-etch step, however, neither adhesive produced shear bond strengths comparable to Clearfil SE, which also produced the most mixed fractures. Storage time did not affect shear bond strengths of any of the materials tested.

## **References**

1. Buoncore MG. A simple method of increasing the adhesion of acrylic filling materials to enamel surfaces. *J Dent Res* 1955; 34:849-53.
2. Van Meerbeek B, Yoshida Y, Van Landuyt K, et al. In: Summitt JB, Robbins JW, Hilton TJ, Schwartz RS, eds. *Fundamentals of Operative Dentistry. A Contemporary Approach*. 3rd edn. Chicago: Quintessence Publishing, 2006:183–260.
3. Nakaayashi N, Pashley DH. *Hybridization of Dental Hard Tissues*. Chicago: Quintessence Publishing; 1998; 65-67.
4. Pereira GD, Paulillo LA, de Ges MF, Dias CT. How wet should dentin be? Comparison of methods to remove excess water during moist bonding. *J Adhes Dent* 2001;3:257-64.
5. Carrilho MRO, Geraldeli S, Tay RF, de Goes MF, Carvalho RM, Tjäderhane L, Reis A, Heling J, Mazzoni A, Breschi L, Pashley DH. In vivo preservation of the hybrid layer by Chlorhexidine. *J Dent Res* 2007; 86:529-33.

6. Van Meerbeek B, De Munck J, Yoshida Y, Inoue S, Vargas M, Vijay P, et al. Buonocore memorial lecture. Adhesion to enamel and dentin: current status and future challenges. *Oper Dent* 2003; 28:215-35.
7. Van Meerbeek B, Yoshihara K, Yoshida Y, Mine A, De Munck J and Van Landuyt KL. State of the art of self-etch adhesives. *Dent Mater* 2010; 1756-68.
8. Sarr M, Mine A, De Munck J, Cardoso MV, Kane AW, Vreven J, Van Meerbeek B, Van Landuyt KL. Immediate bonding effectiveness of contemporary composite cements to dentin. *Clin Oral Investig*. 2010; 14(5):569-77.
9. Heintze SD, Thunpithayakul C, Armstrong ST and Rousson V. Correlation between microtensile bond strength data and clinical outcome of Class V restorations. *Dent Mater* 2011;27:114-125.
10. Peumans M, De Munck J, Van Landuyt K, Lambrechts P, Van Meerbeek B. Five-year clinical effectiveness of a two-step self-etching adhesive. *J Adhes Dent* 2007 Feb;9(1):7-10.
11. Ermis RB, Temel UB, Cellik EU, Kam O. Clinical performance of a two-step self-etch adhesive with additional enamel etching in Class III cavities. *Oper Dent* 2010 Mar-Apr;35(2):147-55.
12. Peumans M, Kanumilli P, De Munck J, Van Landuyt K, Lambrechts P, Van Meerbeek B. Clinical effectiveness of contemporary adhesives: A systematic review of current clinical trials. *Dent Mater* 2005;21:864-881.
13. De Munck J, Van Landuyt K, Peumans M, Poitevin A, Lambrechts P, Braem M, Van Meerbeek B. A critical review of the durability of adhesion to tooth tissue: methods and results. *J Dent Res* 2005;84(2):118-32.
14. Gwinnett AJ, Matsui A. A study of enamel adhesives. The physical relationship between enamel and adhesives. *Arch Oral Biol* 1967;12:1615-9.
15. De Munck J, Van Meerbeek B, Yoshida Y, Inoue S, Vargas M, Suzuki K, Lambrechts P, Vanherle G. Four-year water degradation of total-etch adhesives bonded to dentin. *J Dent Res* 2003;82(2):136-140.
16. Inoue S, Van Meerbeek B, Vargas M, Yoshida Y, Lambrechts P, Vanherle G. Adhesion mechanism of self-etch adhesives. In: Tagami , Toledano M, Prati C. *Advanced adhesive dentistry*. Canada International Symposium 1999. Cirimido (Como): Grafiche Erredue, 2000: 131-148.
17. Sundfeld RH, de Oliveira CH, da Silva AM, Briso AL, Sundfeld ML. Resin tag length of one-step and self-etching adhesives bonded to unground enamel. *Bull Tokyo Dent Coll* 2005;46(3):43-9.
18. Erikson, RL, Barkmeier WW, Kimmes NS. Bond strength of self-etch adhesives to pre-etched enamel. *Dent Mater* 2009;25:1187-1194.

19. Tay FR, Pashley DH, King NM, Carvalho RM, Tsai J, Lai SC, Marquezini L Jr.. Aggressiveness of self-etch adhesives on unground enamel. *Oper Dent* 2004;29(3):309-16.
20. De Munck J, Vargas M, Iracki J, Van Landuyt K, Poitevin A, Lambrechts P, Van Meerbeek B. One day bonding effectiveness of new self-etch adhesives to bur-cut enamel and dentin. *Oper Dent* 2005;30(1):39-49.
21. Van Landuyt KL, Snauwaert J, De Munck J, Peumans M, Yoshida Y, Poitevin A, Coutinho E, Suzuki K, Lambrechts P, Van Meerbeek B. Systematic review of the chemical composition of contemporary dental adhesives. *Biomaterials* 2007;28(26):3757-3785.
22. Hashimoto M, Ohno H, Sano H, Kaga M, Oguchi H. In vitro degradation of resin-dentin bonds analyzed by microtensile bond test, scanning and transmission electron microscopy. *Biomaterials* 2003;24:3795-3803.
23. Hashimoto M, Tay FR, Ohno H, Sano H, Kaga M, Yiu C, Kumagai H, Kudou Y, Kubota M, Oguchi H. SEM and TEM analysis of water degradation of human dentinal collagen. *J Biomed Mater Res B Appl Biomater* 2003;66:287-298.
24. Al-Salehi SK, Burke FJ. Methods used in dentin bonding tests: An analysis of 50 investigations on bond strength. *Quint Int.* 1997;28:717–723.